

A man with a beard and a blue cap is shown in profile, looking towards the left. He is wearing a blue denim jacket and holding a tablet computer. The background is a blurred outdoor setting with trees and a body of water. The text is overlaid on a white diagonal shape.

— THE ORIGINAL —
LASERCAP®

**Prescription-Strength Laser Therapy
For Hair Regrowth**

Regenerate. Regrow. Restore.

Restore your natural hair with the regenerative power of light.

www.LaserCapMD.com



WHAT IS LASERCAP?

LaserCap is the **original** photomedicine cap for hair regrowth. A light and flexible device that can be worn under your favorite hat, LaserCap contains true laser diodes, which illuminate the hair follicles with 650 nm red light. LaserCap is **FDA-cleared** to promote hair regrowth in both men and women with androgenic alopecia.

- **Effective**
Prescription strength laser therapy, FDA-cleared for hair regrowth in men and women.
- **Safe**
Natural, drug-free, and not associated with any serious adverse side effects.
- **Convenient**
Portable, hands-free, and discrete, use LaserCap less than two hours per week at home or on the go.

FDA
CLEARED



INDICATIONS FOR USE


The LaserCap is a 510k FDA-cleared, prescription-only, home use therapeutic device indicated to promote hair growth in men and women with Androgenic Alopecia, Fitzpatrick skin phototypes I to IV, Ludwig-Savin classifications of I-II in women, and Hamilton-Norwood classifications of IIa-V in men.

CONTRAINDICATIONS FOR USE

Patients should not use LaserCap if they have been diagnosed with skin cancer of the scalp, head, or neck.


How does it work?

LaserCap uses **Low-Level-Laser Therapy (LLLT)**, to re-energize inactive hair follicles. Red laser light is absorbed by specialized cells deep within the follicle, providing energy to enhance cellular respiration and hair growth. As a result, follicles become awakened from their inactive state and produce thicker, healthier hair. LLLT is a safe and effective treatment for both male and female pattern hair loss, proven in numerous high quality clinical trials.

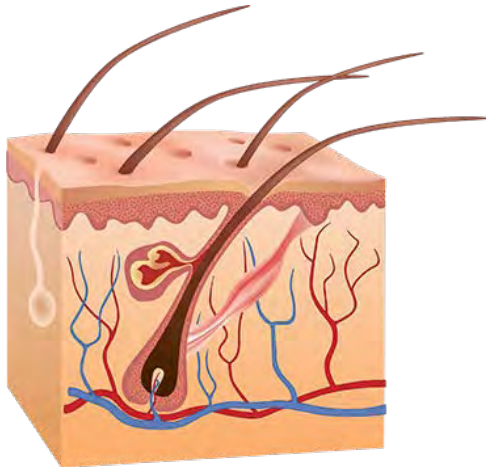


Before

Follicles are inactive, growing slowly and producing thin hair.

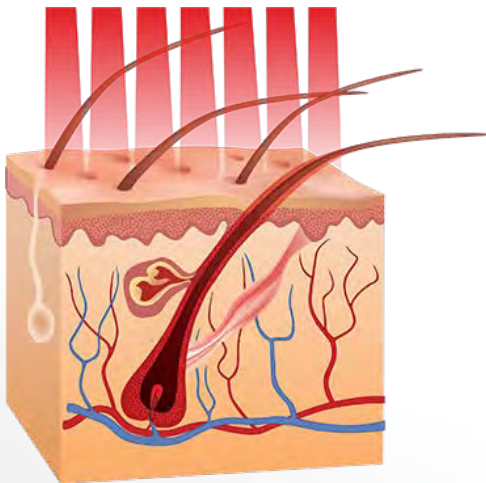

During

Laser light is absorbed by the cells of the hair follicle, providing energy to enhance cellular respiration and hair growth.

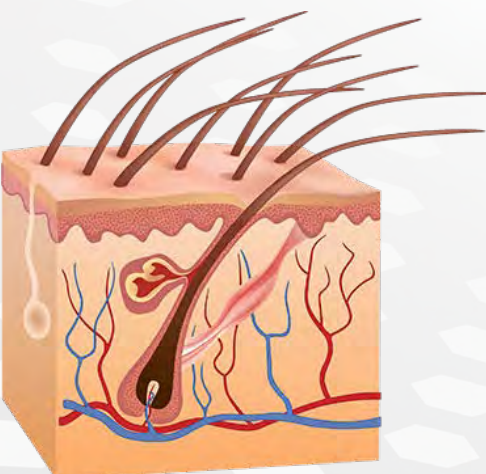

After

Dormant follicles are awakened from their inactive state and produce thicker, healthier hair.


PRESCRIPTION STRENGTH LASER THERAPY FOR HAIR REGROWTH



Did You Know?

LaserCap delivers a **higher dose of energy** per treatment than other comparable devices. Many of these devices in fact deliver a dose of energy well below what experts consider to be optimal.

DOSE OF ENERGY PER TREATMENT (J/cm²)

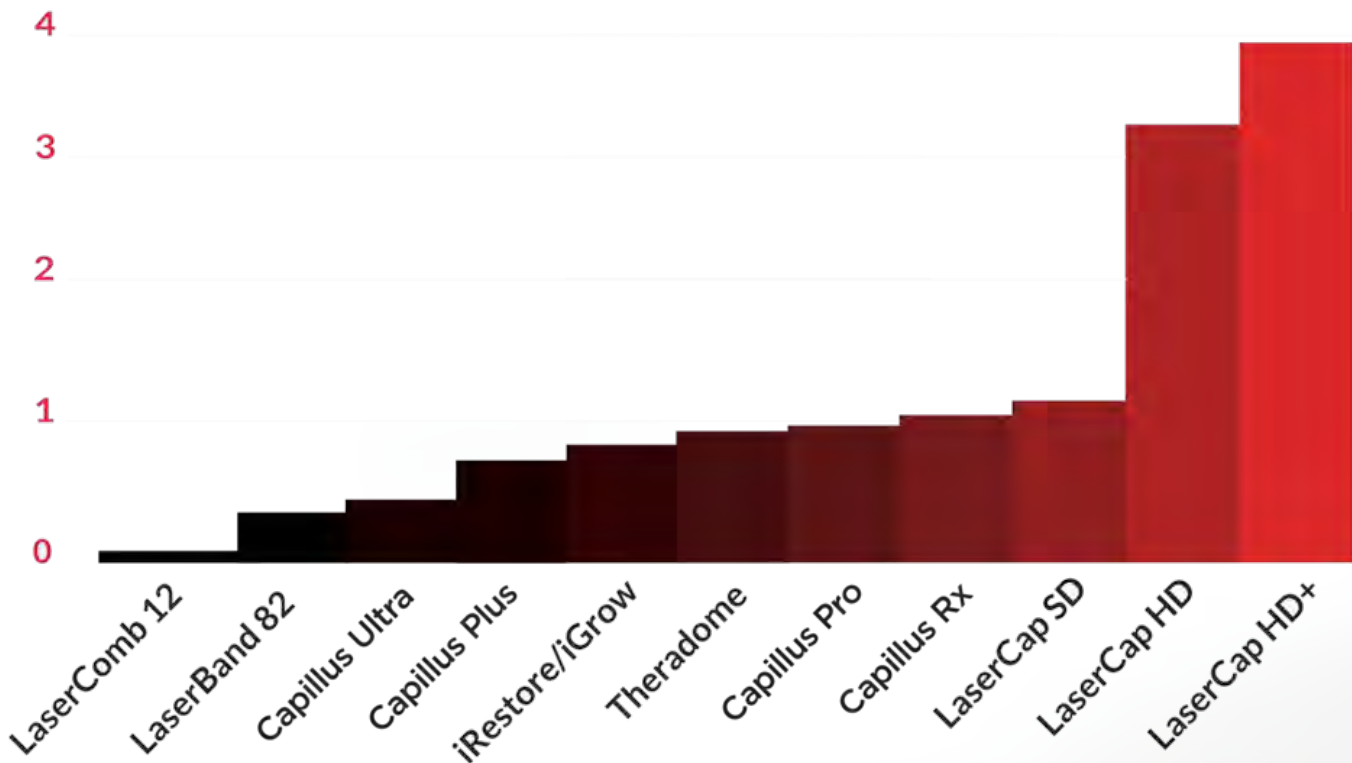

Pulsed Wave Emission

Key to LaserCap's prescription-strength treatment is **pulsed wave light emission**. This is a fancy way of saying LaserCap turns off and on its lasers many times per second during a treatment, causing them to blink, or flicker. Doing so allows the laser diodes to cool during the fraction of a second when they switch off. As a result, LaserCap is able to deliver higher levels of energy without excessive build-up of heat.

In addition to **prescription-strength**, LaserCap is uniquely designed to maximize **ease of use** and **treatment consistency**:

- **PRESCRIPTION-ONLY**
Prescription-strength laser therapy monitored by a medical professional.
- **FULL COVERAGE**
The entire affected area of the scalp is illuminated at the same time, improving consistency of treatments.
- **HANDS-FREE**
Rest and relax during treatments instead of combing through your hair or otherwise moving the device around.
- **PORTABLE**
Powered by a rechargeable battery pack, use LaserCap at home or on the go.
- **DISCREET**
Not a big, bulky helmet, but instead can be worn discreetly under your favorite hat.

| LASERCAP | THERADOME | iRESTORE |
|---------------------|---------------------|---------------------|
| ✓ Prescription-Only | ✗ Prescription-Only | ✗ Prescription-Only |
| ✓ Full Coverage | ✓ Full Coverage | ✓ Full Coverage |
| ✓ Hands-Free | ✓ Hands-Free | ✓ Hands-Free |
| ✓ Portable | ✓ Portable | ✓ Portable |
| ✓ Discreet | ✗ Discreet | ✗ Discreet |

| iGROW | LASERBAND | LASERCOMB |
|---------------------|---------------------|---------------------|
| ✗ Prescription-Only | ✗ Prescription-Only | ✗ Prescription-Only |
| ✓ Full Coverage | ✗ Full Coverage | ✗ Full Coverage |
| ✓ Hands-Free | ✗ Hands-Free | ✗ Hands-Free |
| ✗ Portable | ✓ Portable | ✓ Portable |
| ✗ Discreet | ✗ Discreet | ✗ Discreet |

A LIFETIME OF VALUE

The LaserCap is built to last.

Did you know we have patients who have used the same LaserCap for over 10 years?

LaserCap Company is so confident in the quality of our products, we are proud to offer a best-in-industry **lifetime warranty policy** and **one-year satisfaction guarantee**. With our world-class customer service, you can be assured a lifetime of therapy to regrow and restore your natural hair.




The LaserCap Lifetime Warranty Policy covers defects in materials or workmanship. If a LaserCap is determined to be defective within 3 years from date of purchase, LaserCap company will repair or replace for a like model at no charge. After 3 years from date of purchase a nominal replacement fee will be charged. The one year satisfaction guarantee allows LaserCap to be returned within 1 year from date of purchase minus a 25% return fee and along with shipping charges, credit card fees, taxes, and any additional fees incurred by LaserCap Company or their authorized physician partners.




Real People. Real Results.

LaserCap has helped thousands of patients regrow and restore their natural hair. The following photos show just a few examples of the results that can be obtained using LaserCap.

These are all **real** LaserCap patients, not models, from case studies conducted by LaserCap physicians.

All follow-up photos were taken after 4-6 months of treatment with LaserCap only.

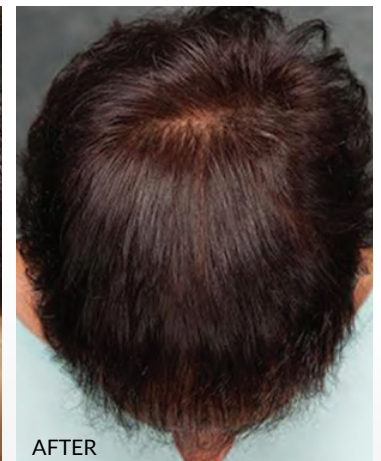

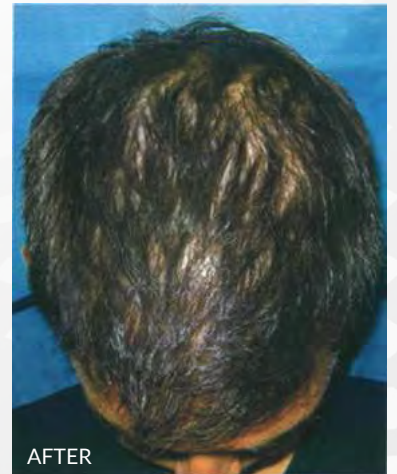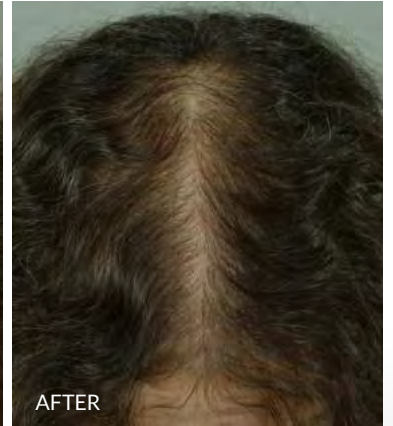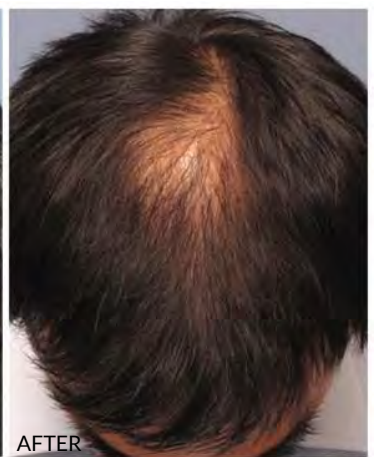

FAQs

How often do I use LaserCap?

Our recommended treatment protocol is to wear the cap for 30 minutes every other day. Your provider however may recommend an alternative protocol.

How long before I see results?

If you are compliant with your treatment, typically you will have noticeable improvement in 4 to 6 months.

Are there any side effects?

LaserCap is not associated with any serious adverse side effects. A small percentage of LaserCap users have reported a mild headache, and minor pruritus of the scalp (an itchy or tingling sensation). If experienced at all, these side effects are always transient, occurring only during use.

Can both men and women use LaserCap?

Yes! LaserCap is equally effective in both men and women.

Can I use LaserCap with other treatments for hair loss?

Yes! LaserCap is effective in combination with all other treatments for hair loss. Your provider may recommend a combination therapy for improved results.



Warning!



Is that "laser cap" a lemon? A warning about cheap LaserCap copies.

Did you know that The Original LaserCap is one of the most widely copied technologies in the world? Many of these copy devices can be found online for low prices. While some may look very similar to LaserCap, they should be avoided at all costs, for several reasons:

- **Low-Strength**
LaserCap copies will likely not be effective because they have not optimized their output, and deliver light energy at very low intensity.
- **Poor Quality**
Copy devices are often made 100% in China with poor quality control. Frequently they malfunction or break down after a short period of use.
- **Limited Warranties**
No LaserCap copies have a lifetime warranty policy. Often they have very limited warranty, or no warranty at all.

Nothing beats the original.

Don't waste your money on cheap "laser cap" lemons. Stick with **The Original LaserCap** instead!

- **FDA-cleared technology**
- **Prescription-strength treatment**
- **Backed by experts**
- **Lifetime warranty policy**
- **Assembled in USA**



— THE ORIGINAL —
LASERCAP[®]

Prescription-Strength Laser Therapy
For Hair Regrowth

CONTACT US

www.LaserCapMD.com

info@lasercapmd.com

1-855-424-7774

LaserCap Company
26 Alpha Park Dr.
Highland Heights, OH
44143